


How To Assemble Your DIY Face Shield V11

Paper instructions below. [Video here.](#)

This face shield is cleanable, reusable. The shield can be reused, soaked in Isopropyl Alcohol/Isopropanol 70% for 5 minutes or Sodium Hypochlorite (household bleach) 0.01% (1:10), 5 mins+.¹

Parts included in your kit

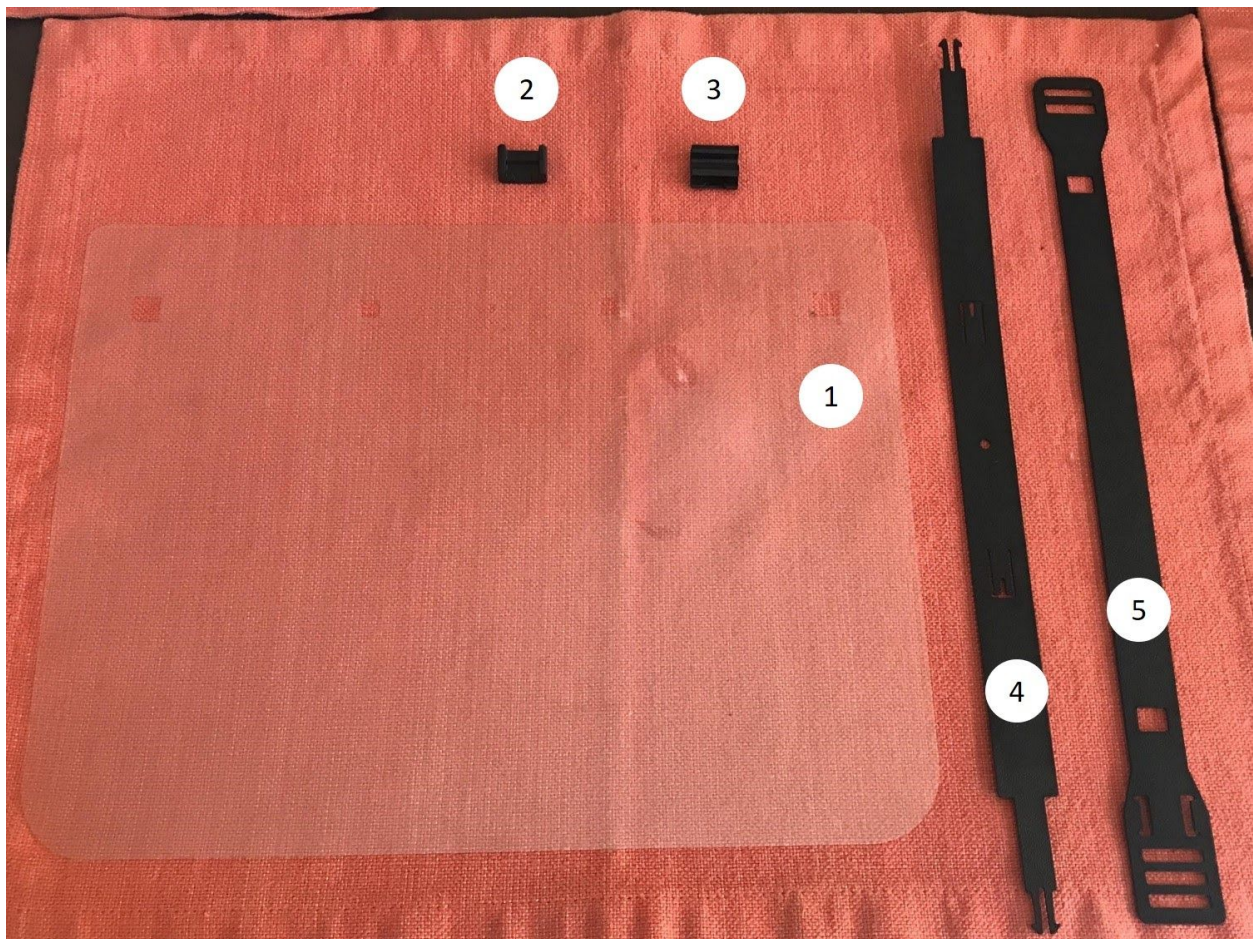

¹

https://help.prusa3d.com/en/article/prusa-face-shield-disinfection_125457?fbclid=IwAR2DXQi11XWHamiwdIKqQanED50CPu4jKb40f-3TesGO-3XLkPiSG0ygDdQ

#1 - Shield

#2 - 1 Small Clip

#3 - 1 Large Clip


#4 - Pronged band

#5 - Slotted band

Steps for complete assembly

Note: Your kit may have arrived partially assembled. If so, remove the clip attached to the 2-slotted end of the slotted band (#5) by pushing on the ends of the clip so that it pops out. Then skip to Step 6.

Step 1 The shield has 4 holes. Weave the pronged band (#4) into the 4 holes of the shield, first the outer holes, and then wiggle the black pieces into the inner holes. This band will be on the outside while the prongs will protrude inside.

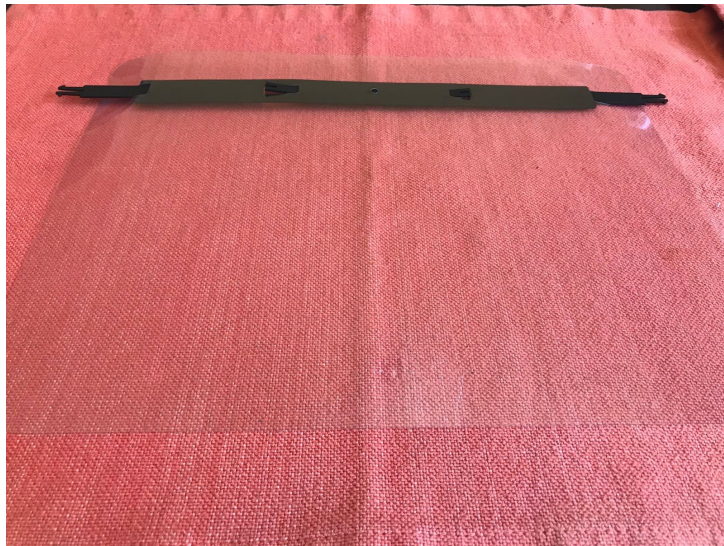


Step 2 With the inside of the shield (prong side) facing you, put one of the pronged ends through the square hole on the 3-slotted end of the slotted band (#5).



Step 3 Take the thinner of the 2 clips (#2) and attach it beneath the prongs on the 3-slotted end, with the groove of the clip facing the shield and the back side facing you.



Step 4 Insert the prongs into the clip (#2) until you hear a click. To release, squeeze the prongs and pull out. Note: These materials are brittle and disassembly is not recommended.

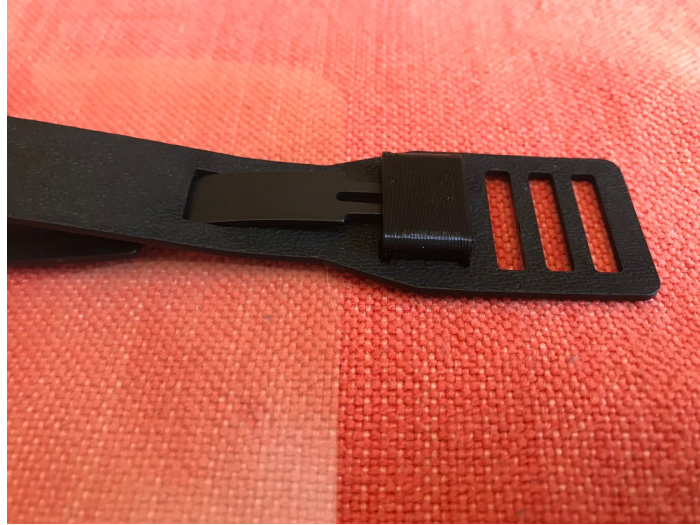

Step 5 Weave the remaining pronged end through the square hole on the 2-slotted end of the slotted band (#5).



Step 6 On the 2-slotted end, snake the shoestring (or elastic) through the 2 slots, starting and ending at the front of the shield. Take the remaining clip (#3) and attach it over the slots and shoestring, with the groove facing the shield and the back side facing you. Push until you hear a click.



Step 7 Insert the prongs into the clip (#3) and push until you hear a click.



Step 8 Snake the shoestring into the first slot on the 3-slotted side of the band, then go out the third slot, and return into the second slot, finally going back out again through the first slot. The length may be adjusted by pulling the string. Now your face shield is complete!





These materials are courtesy of [DIYPPE.org](https://www.diypppe.org)

Our thanks go out to the healthcare workers for being on the front lines of this crisis. You are our relatives, friends and neighbors.

These materials and our website are for you, our heroes.

IMPORTANT NOTICE: These Face Shields have not been tested or approved by any governmental agency and are not offered with any promise of safety or efficacy in protecting against any disease or infection, including protection against COVID-19. Face Shields ARE NOT PROVIDED FOR ANY PARTICULAR PURPOSE. THE PROVIDER AND ALL SUPPLIERS MAKE the Face Shields AVAILABLE free of charge AS-IS AND WITH ALL FAULTS AND DISCLAIM ALL WARRANTIES, EXPRESS OR IMPLIED, AS TO QUALITY OR EFFECTIVENESS FOR ANY PURPOSE. THE USER ASSUMES ALL RISK OF USE.

WARNING: Some products may contain chemicals known to the State of California to cause cancer, or birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov.